



CONTACT US

LOCATIONS

561 E Garden Drive, Unit H
Windsor, Colorado 80550
970-686-2266

1750 N. Humboldt St. #102
Denver, Colorado 80218
303-997-6753

2110 Evans Ave
Cheyenne, WY 82001
307-778-6595

*TAKE CONTROL OF YOUR FIT
WITH THE QUATRO!*

QUATRO LANYARD

KISS LANYARD
CLICK LANYARD

QUATRO PIN-LOCK

BULLDOG APL
OSSUR ICE LOCK
FILLAUER ORIGINAL PIN-LOCK
BULLDOG DROP-IN
OWW ALPHA LOCK

QUATRO SUCTION

LYNN VALVE SMALL
LYNN VALVE LARGE

QUATRO

ADJUSTABLE
SOCKET

Quorum
3D LABS & PROSTHETICS



3D SCAN

All 3D printed sockets start with a scan of the residual limb, a plaster cast or a well fitting socket.





ENGINEER

Our highly skilled team of engineers will use the digital scan to make the socket 100% custom for the individual.




PRINT


Each print is additively manufactured using the highest quality materials and state of the art technology.

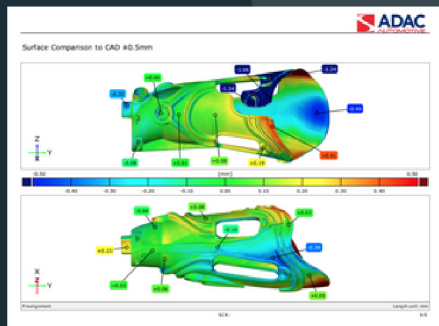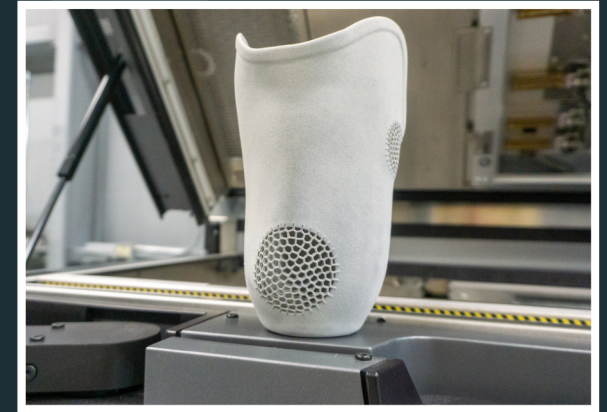

Quatro sockets are rigorously tested to ensure you will have a lightweight and incredibly strong socket for many years.

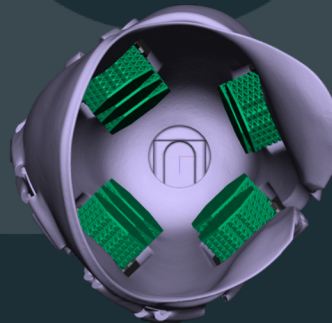

Quatro sockets are designed with four zones of compression control. Each zone offers up to 60lbs of grip force and up to 12% volume adjustments using an adjustable dial.

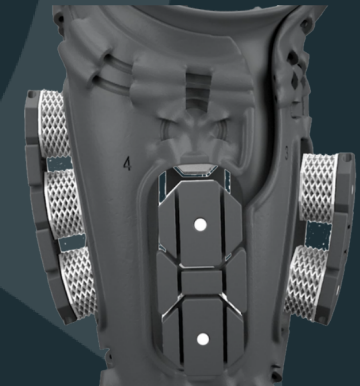

All of these components work in unison to improve proprioception, resulting in less falls and more freedom.